


Updated 2019

**hpi**

Hypoxyprobe, Inc.  
121 Middlesex Turnpike Burlington,  
MA 01803 USA  
[www.hypoxyprobe.com](http://www.hypoxyprobe.com)

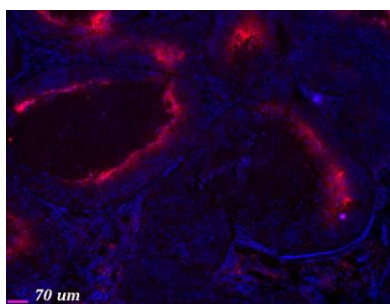

## Data Sheet: Hypoxyprobe™ RedAPC Kit (Cat# HP8-x)

Description:	Mouse IgG1 monoclonal antibody conjugated to Allophycocyanin (APC) fluorophore (HP-RedAPC-MAb) and x mg of solid pimonidazole HCl (x = 100, 200 or 1000 mg).
Specificity:	Pimonidazole forms adducts with thiol containing proteins at $pO_2 \leq 10$ mm Hg. HPRedAPC-MAb binds to pimonidazole adducts in hypoxic cells in tissues and in culture. Please see <a href="http://www.hypoxyprobe.com">www.hypoxyprobe.com</a> for details of the Hypoxyprobe system for detecting cell and tissue hypoxia.
Format:	Each vial of HP-RedAPC-MAb contains anti-pimonidazole antibody dissolved in 200 microliters of stabilized buffer at a concentration of approximately 0.50 mg/mL.
Product Type:	Fluorophore conjugated IgG <sub>1</sub> mouse monoclonal antibody (clone 4.3.11.3)
Fluorophore-Protein Ratio:	$\geq 1$ APC moiety per IgG <sub>1</sub> molecule.
Fluorescence:	Excitation 633 nm; Emission max 660 nm. Allophycocyanin has strong red fluorescence.

## Product Details

Applications:	Flow cytometry on isolated cells and immunofluorescence on frozen tissue sections. HP-APC-MAb is an alternative to FITC-labeled anti-pimonidazole mouse monoclonal antibody <sup>1</sup> . 1/20-1/200 dilution is suggested but users should optimize in their system using appropriate negative and positive controls. Image above: human tumor xenograft hypoxia (red); 1/200 dilution of HP-RedAPC overnight at 4°C on frozen section. (Courtesy of Hans Peters, Radboud University Nijmegen Medical Centre, Nijmegen, The Netherlands.)
Target Species:	All.
Product Form:	Purified IgG <sub>1</sub> prepared by affinity chromatography
Buffer:	Proprietary formulation containing stabilizers.

1. Jankovic B, Aquino-Parsons C, Raleigh JA, et al. Comparison between pimonidazole binding, oxygen electrode measurements, and expression of endogenous hypoxia markers in cancer of the uterine cervix. *Cytometry B Clin Cytom* 2006; 70: 45-55.